

EASY ADJUST MOBILE LEAD PORTA SHIELD

PRODUCT INFORMATION

PHILLIPS-SAFETY.COM

COPYRIGHT ©2026. ALL RIGHTS RESERVED



PRODUCT INFORMATION

EASY ADJUST MOBILE LEAD PORTA SHIELD



The Easy Adjust Mobile Lead Porta-Shield is the ultimate solution for clinicians who require effortless, highly adjustable radiation protection. This model replaces traditional manual-lift shields with a high-torque hand crank system, designed to take the physical strain out of height adjustments. Simply turn the handle to glide the shield up or down—no heavy lifting required.

Key Features:

- 1. Effortless “Stop-Anywhere” Control:** The fluid-motion hand crank allows you to dial in the exact height you need. The internal mechanism holds the shield securely in place the moment you stop turning, providing infinite adjustability without the hassle of pins or preset holes.
- 2. Ergonomic Design:** Engineered for ease of use, the mechanical drive system does the heavy lifting for you, allowing for smooth travel between 24 and 48 inches with minimal physical effort.
- 3. Superior Mobility:** Mounted on a heavy-duty five-caster base, the 24-inch diameter footprint ensures a sturdy, tip-resistant foundation that glides effortlessly across hospital floors.
- 4. Certified 0.50 mm Pb Protection:** Delivers industry-standard medical X-ray attenuation, offering a simple yet effective barrier for personnel and patients alike.

SPECIFICATION

COLOR: Black/Blue


LEAD PROTECTION: 0.50 mm Pb Eq

SIZE: 24” x 24” | 24” x 36” | 24” x 48”

FEATURES:

Speaker-Crank Drive, Lock-Tight Geometry,
All-Terrain Casters, Lead Core

PRODUCT INFORMATION
EASY ADJUST MOBILE LEAD PORTA SHIELD

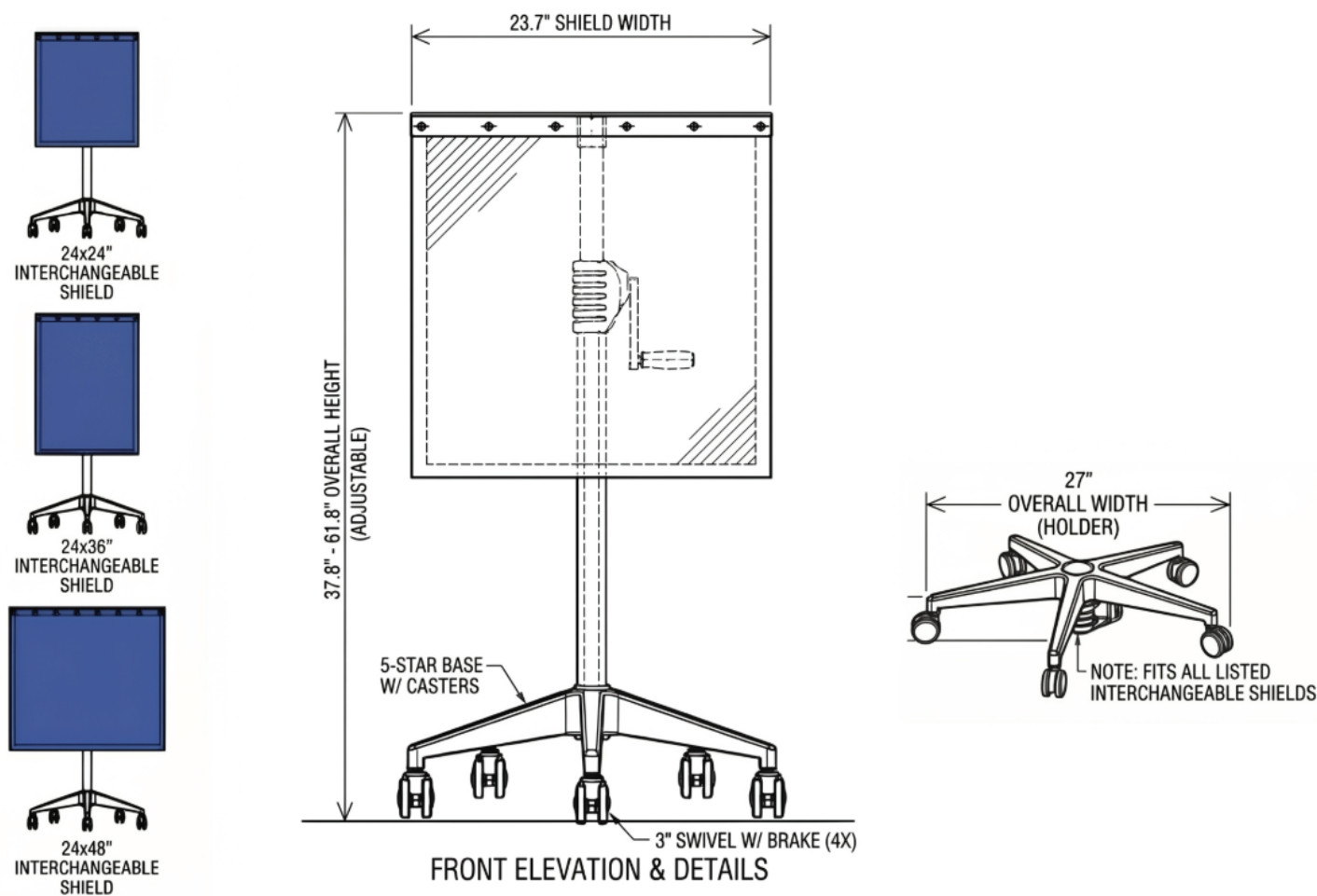

Date: 6/24/2026		Drawing Title: Easy Adjust Mobile Lead Porta Shield
Scale: 1:10	Unit: Inches	
		Drawing No: LB-MB2424-HC, LB-MB2436-HC, LB-MB2448-HC
		Notes: All dimensions $\pm 0.1"$

CONTACT

Should you need any further information,
please do not hesitate to contact us.

123 Lincoln Blvd, Middlesex, NJ 08846, USA 

+1 866 575 1307 

service@phillips-safety.com 

www.phillips-safety.com 